



# Kingon P2-E

## Portable Oxygen Concentrator

Reliable pulse dose oxygen therapy designed for mobility without compromising performance.



### Lightweight & Portable

Approx. 1.98 kg compact design for everyday mobility.



### Intelligent Pulse Dose Delivery

Breath-triggered oxygen support with settings 1–5



### Quiet Operation

Low noise output suitable for home, travel and public settings.



### Extended Battery Operation

Up to 5 hours per charge, depending on setting.




### Travel Ready Power System

FAA approved with rechargeable battery and AC/DC charging

# Designed for Reliable Supplemental Oxygen Therapy On The Go

The Kingon P2-E Portable Oxygen Concentrator provides pulse dose oxygen therapy in a compact, lightweight format suited for individuals wanting portable oxygen.

Engineered for consistent performance, the P2-E supports mobility without compromising delivery reliability.

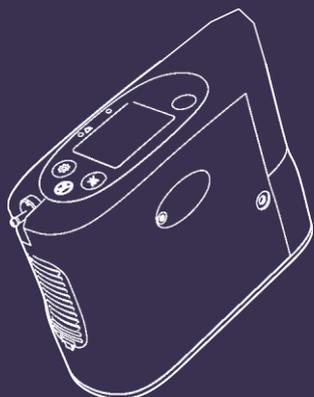

## Intelligent Oxygen Delivery & Performance

The P2-E delivers breath-triggered pulse dose oxygen while maintaining stable concentration across all settings.

Pulse dose flow settings: 1 - 5 Intelligent inhalation detection Smart auto-pulse mode

### Performance Output

Oxygen concentration: 90% -3%/+6%  
Stable output across pulse settings  
Low noise operation suitable for daily use



## Impressive Battery Duration

The Kingon P2-E features a rechargeable lithium battery designed to support mobile oxygen therapy across various pulse settings.

Pulse Setting	Battery Duration
Flow Setting 1	5 hours
Flow Setting 2	4 hours
Flow Setting 3	3.25 hours
Flow Setting 4	2.5 hours
Flow Setting 5	1.7 hours



### Integrated Safety monitoring

Built-in alarms provide operational reassurance, alerting you to low oxygen concentration, low battery, no breath detected, or system faults.



### Travel-Ready Power System

FAA approved for air travel, the Kingon P2-E includes an AC adaptor, DC car charger, and a rechargeable lithium battery



# Built for Everyday Mobility

Weighing approximately 1.98 kg, the Kingon P2-E is designed to support active use across home, community and travel environments.

Its compact footprint and included carry bag allow for convenient daily transport while maintaining reliable oxygen delivery.

## Compact & Lightweight Design

Engineered for portability, the Kingon P2-E supports mobile oxygen therapy in a compact, travel-friendly format.


Weight: Approx. 1.98 kg  
Dimensions: 22 cm (H) x 16 cm (W) x 8.5 cm (D) Ergonomic carry design with included carry bag

Its lightweight construction supports comfortable extended use across home, community and travel environments.



## Clear, Simple Operation

The clear LCD display and straightforward controls make the unit easy to use, even for first-time oxygen concentrator users.

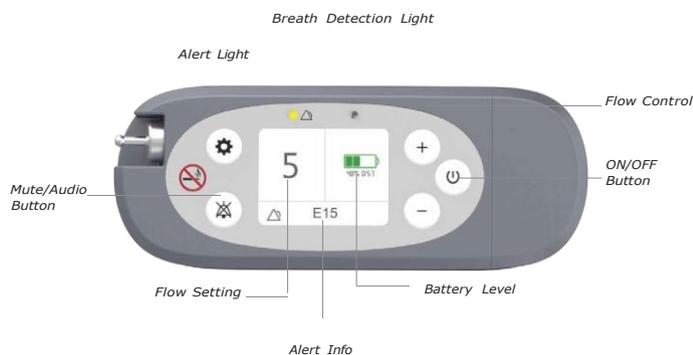


### Low Noise Operation

Kingon P2-E operates at just 37 dB(A), keeping use discreet in shared or public spaces



### 5 Years Total Warranty

Your Kingon Portable Oxygen Concentrator includes a 3-year warranty, plus a 2-year warranty upgrade when registered on [kingon.com.au](http://kingon.com.au)

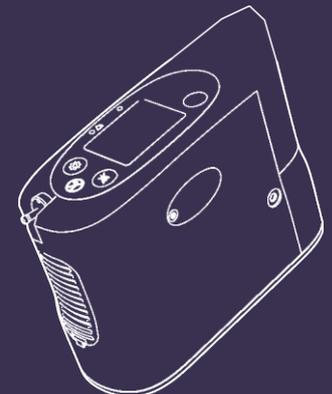

# Technical Specifications

Weight	1.98kg (includes battery)
<b>Size</b>	Length: 221mm
	Width: 85mm
	Height: 160mm (includes battery)
<b>ARTG</b>	Included in the ARTG

## Flow Setting & Pulse Volume

Oxygen Flow	Pulse dose delivery system
Oxygen Concentration	90% -3%/+6% at all settings
Warm-Up Time	2 minutes
Display Screen	2.8 inch

## Battery & Power

Battery	Five flow settings: 1, 2, 3, 4, 5
	Duration: up to 5 hours
	Recharge: approximately 4 hours with AC/DC power supply

<b>Power</b>	AC power supply 100-240V, 50-60Hz (Auto sensing to allow worldwide use)
<b>Voltage:</b>	14.4VDC, rated capacity: 6.8Ah
Sound Level	37 dB(A) (on setting 2) (Prüfmethode 141 10/2018 MDS-Hi)

## Warranty & Quality Standards

Quality Standards	FDA, CE, FAA Approved
In Flight Use	Approved
Warranty	(3-year standard + 2-year free extension when you register your device) 1 year for battery & sieve bed
<b>Operation</b>	Operation Simple control functions and easy-to-read LCD display

## Package Includes:



Distributor:



Toll Free: 1800 812 097  
E: [sales@wymedical.com.au](mailto:sales@wymedical.com.au) W: [wymedical.com.au](http://wymedical.com.au)  
466 Fullarton Road, Myrtle Bank SA 5064